


V-Curette sterile, disposable

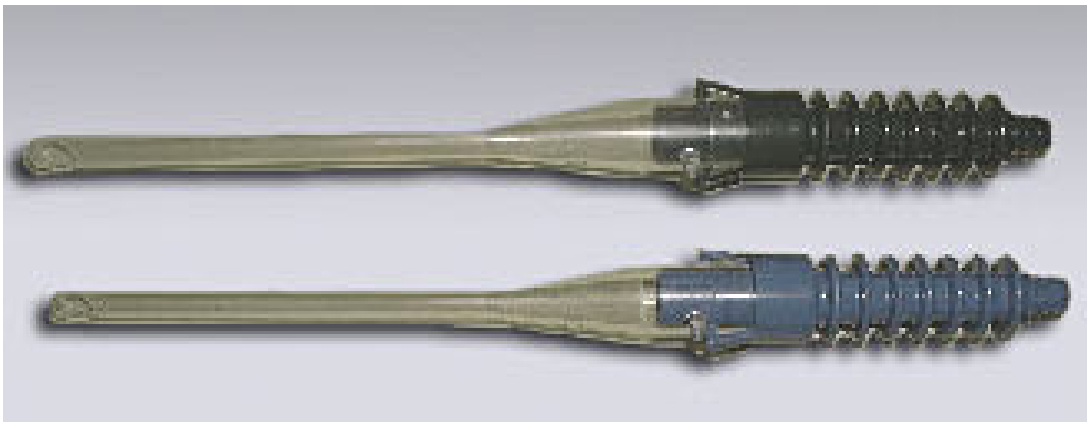

The Curette and the 1,5 metre plastic tubing with tube connector are sterilized by electron beam radiation.

The V-Curette is available in two sizes. Blue handle 8 mm outside diameter, green handle 10 mm. The tubing used for both sizes is identical.

REF 24080: V-Curette 8 mm diameter, complete with handle.

REF 24100: Ditto, 10 mm diameter.

REF 24110: 1,5 m PVC V-Tube with tube connector.

The V-Curette made of transparent plastic (Cyrolite) can be rotated to a desired position by means of the knobs in front of the handle. With the proper technique it is very often possible to use the 8 mm Curette up to twelve weeks of pregnancy, thus eliminating unnecessary dilation and minimizing the blood loss.

Due to the design of the suction slot, the operator can feel the surface of the uterine cavity so well that there is no need of other instruments to make sure that the cavity has been completely emptied.